

Document information

Software: realplex 2.2

File Name: EPPENDORF\Lorenzo\QPCR16.04.2

Printed by: EPPENDORF
Created: Apr/16/2018 12:39

Serial No. Thermo Module: 6325 30387 Serial No. realplex Module.: 630011465

Acquisition Start Time: EPPENDORF Apr/16/2018 12:45
Acquisition End Time: EPPENDORF Apr/16/2018 14:13
Last updated: EPPENDORF Apr/06/2018 12:28

Background: Sarstedt-20µl Sep/12/2011 10:28 Color Calibration: SYBR Mar/12/2018 15:31

QPCR16.04.2018 Quantification Apr/16/2018 14:13

Inverted Data: OFF

Comment:

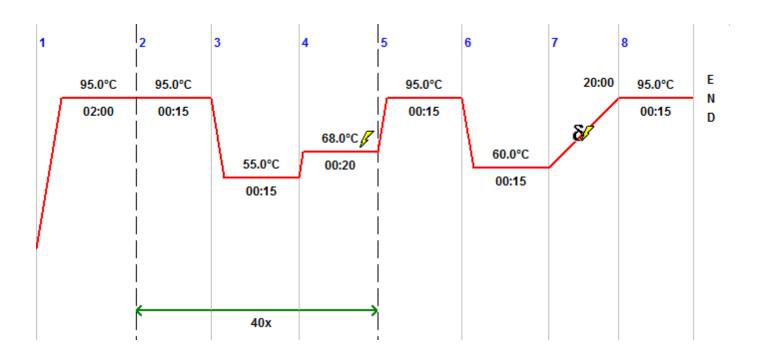


Plate layout

	1	2	3	4	5	6	7	8	9	10	11	12
A	CEWE	CEWE	CEWE	CEWE	CEWE	CEWE	ILWE_A	ILWE_A	ILWE_A	ILWE_A	ILWE_A	ILWE_A
	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00
В	CEWE	CEWE	CEWE	CEWE	CEWE	CEWE	ILWE_A	ILWE_A	ILWE_A	ILWE_A	ILWE_A	ILWE_A
	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00
С	CEWE	CEWE	CEWE	CEWE	CEWE	CEWE	ILWE_A	ILWE_A	ILWE_A	ILWE_A	ILWE_A	ILWE_A
	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00
D	CEWE	CEWE	CEWE	CEWE	CEWE	CEWE	ILWE_A	ILWE_A	ILWE_A	ILWE_A	ILWE_A	ILWE_A
	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00
E	CEWE	CEWE	CEWE	CEWE	CEWE	CEWE	ILWE_A	ILWE_A	ILWE_A	ILWE_A	ILWE_A	ILWE_A
	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00	1: 1.00
F												
G												
Н	NTC	NTC	NTC	NTC	NTC	NTC	water	water	water	water	water	water



PCR Program




Program Header

Lid Temp	105 °C	TSP Heated Lid	Yes
Temp. Mode	Standard	Switch off lid at low block temp	No
Impulse	No	Simulate Mastercycler gradient	No




Raw Data SYBR

Fluorescence Profile

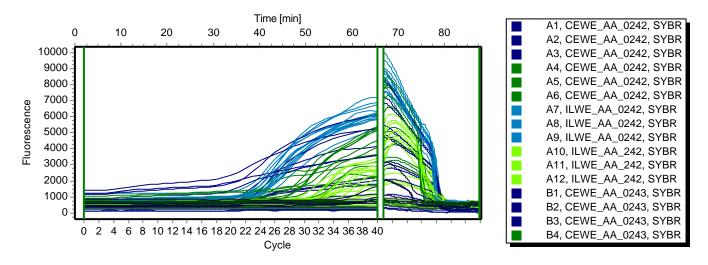

Temperature Profile

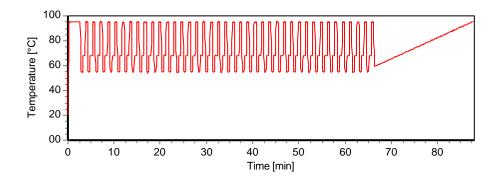



Quantification SYBR


Pos	Name	Ct SYBR	Ct Mean SYBR	Ct Dev. SYBR	Amount SYBR [Copies]	Amount Mean SYBR	Amount Dev. SYBR	Target SYBR
! ■ A1	CEWE_AA_0242	38.49			1.00			mouse
! ■ A2	CEWE_AA_0242				1.00			mouse
! ■ A3	CEWE_AA_0242				1.00			mouse
! ■ A4	CEWE_AA_0242				1.00			eimeria
! ■ A5	CEWE_AA_0242				1.00			eimeria
! ■ A6	CEWE_AA_0242				1.00			eimeria
! □ ►A7	ILWE_AA_0242	21.67	22.52	1.04	1.00			mouse
! ■ A8	ILWE_AA_0242	22.22	22.52	1.04	1.00			mouse
! ■ A9	ILWE_AA_0242	23.68	22.52	1.04	1.00			mouse
! 	ILWE_AA_242				1.00			eimeria
! 	ILWE_AA_242	37.49			1.00			eimeria
! ■ A12	ILWE_AA_242				1.00			eimeria
! ■ B1	CEWE_AA_0243		23.19	0.67	1.00			mouse
! ■ B2	CEWE_AA_0243	23.67	23.19	0.67	1.00			mouse
! ■ B3	CEWE_AA_0243	22.72	23.19	0.67	1.00			mouse
! ■ B4	CEWE_AA_0243	21.93	22.15	0.45	1.00			eimeria
! ■ B5	CEWE_AA_0243	21.85	22.15	0.45	1.00			eimeria
! ■ B6	CEWE_AA_0243	22.66	22.15	0.45	1.00			eimeria
! ■ B7	ILWE_AA_0243	22.31	23.04	0.64	1.00			mouse
! ■ ■88	ILWE_AA_0243	23.38	23.04	0.64	1.00			mouse
! ■ B9	ILWE_AA_0243	23.44	23.04	0.64	1.00			mouse
!	ILWE_AA_0243	32.68	31.96	0.89	1.00			eimeria
!	ILWE_AA_0243	32.24	31.96	0.89	1.00			eimeria
! ■ B12	ILWE_AA_0243	30.97	31.96	0.89	1.00			eimeria
! Ⅲ ■C1	CEWE_AA_0247		24.14	2.58	1.00			mouse
! ■ C2	CEWE_AA_0247	25.96	24.14	2.58	1.00			mouse
!	CEWE_AA_0247	22.32	24.14	2.58	1.00			mouse
! 	CEWE_AA_0247	30.90	30.91	0.53	1.00			eimeria
! ■ C5	CEWE_AA_0247	31.44	30.91	0.53	1.00			eimeria
!	CEWE_AA_0247	30.38	30.91	0.53	1.00			eimeria
! ∏	ILWE_AA_0247	21.88	22.05	0.68	1.00			mouse
!	ILWE_AA_0247	21.47	22.05	0.68	1.00			mouse
i	ILWE_AA_0247	22.80	22.05	0.68	1.00			mouse
! ∏	ILWE_AA_0247	31.72	32.33	1.28	1.00			eimeria
! ∏	ILWE_AA_0247	31.48	32.33	1.28	1.00			eimeria



Pos	Name	Ct SYBR	Ct Mean SYBR	Ct Dev. SYBR	Amount SYBR [Copies]	Amount Mean SYBR	Amount Dev. SYBR	Target SYBR
. C12	ILWE_AA_0247	33.80	32.33	1.28	1.00			eimeria
D1	CEWE_AA_0248	20.04	20.35	1.14	1.00			mouse
D2	CEWE_AA_0248	19.40	20.35	1.14	1.00			mouse
! ■ D3	CEWE_AA_0248	21.62	20.35	1.14	1.00			mouse
! ■ D4	CEWE_AA_0248	26.63	27.23	0.52	1.00			eimeria
! ■ D5	CEWE_AA_0248	27.56	27.23	0.52	1.00			eimeria
! ■ D6	CEWE_AA_0248	27.51	27.23	0.52	1.00			eimeria
! ■ D7	ILWE_AA_0248	23.24	22.46	0.73	1.00			mouse
! ■ D8	ILWE_AA_0248	22.35	22.46	0.73	1.00			mouse
! ■ D9	ILWE_AA_0248	21.80	22.46	0.73	1.00			mouse
! □ D10	ILWE_AA_0248	32.90	33.81	1.15	1.00			eimeria
! ☐ D11	ILWE_AA_0248	33.43	33.81	1.15	1.00			eimeria
! □ D12	ILWE_AA_0248	35.10	33.81	1.15	1.00			eimeria
! ■ E1	CEWE_AA_0249	19.11	18.81	0.51	1.00			mouse
! ■ E2	CEWE_AA_0249	19.09	18.81	0.51	1.00			mouse
! ■ E3	CEWE_AA_0249	18.21	18.81	0.51	1.00			mouse
! ■ E4	CEWE_AA_0249	28.72	28.64	0.83	1.00			eimeria
! ■ E5	CEWE_AA_0249	29.42	28.64	0.83	1.00			eimeria
! ■ E 6	CEWE_AA_0249	27.78	28.64	0.83	1.00			eimeria
! ■ E 7	ILWE_AA_0249	21.00	21.86	0.88	1.00			mouse
! ■ E8	ILWE_AA_0249	21.83	21.86	0.88	1.00			mouse
! ■ E9	ILWE_AA_0249	22.75	21.86	0.88	1.00			mouse
! ■E10	ILWE_AA_0249	30.37	29.91	1.08	1.00			eimeria
! ■E11	ILWE_AA_0249	30.69	29.91	1.08	1.00			eimeria
! ■ E12	ILWE_AA_0249	28.68	29.91	1.08	1.00			eimeria
- □ ■H1	NTC	37.40	38.47	1.52	-			mouse
- □ ■H2	NTC	-	38.47	1.52	-			mouse
- □ ■H3	NTC	39.54	38.47	1.52	-			mouse
- □ ■H4	NTC	-			-			eimeria
- □ ■H5	NTC	-			-			eimeria
- □ ■H6	NTC	-			-			eimeria
H7	water	-			-			mouse
- □ ■H8	water	-			-			mouse
- □ ■H9	water	36.49			-			mouse
-U H10	water	-			-			eimeria
- ■H11	water	-			-			eimeria
- ■ H12	water	-			-			eimeria




Amplification Plot

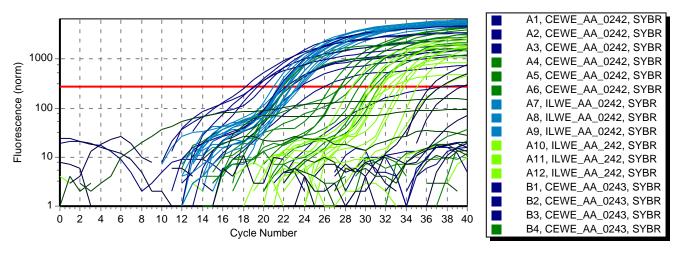

Threshold 276 (Noiseband)

Baseline automatic, Drift correction OFF

Standard curve

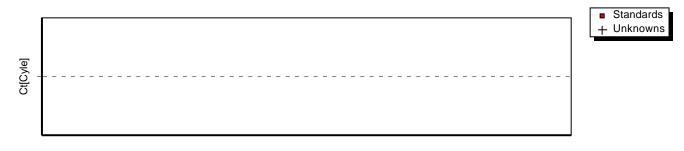

Amount[Copies]

Slope - R^2 -Y-Intercept - Efficiency -